


## INVERSE ULTRASONIC OSMOSIS FOR FAT LOCALIZED ADIPOSITY



Intuitive touch screen with preset protocols and customizable programs. Software available in all languages.

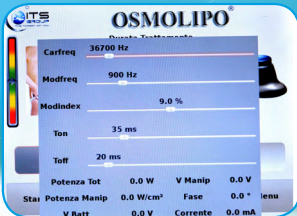


Phase control device to ensure the effectiveness of treatment.



Painless body treatment:

- Fast sessions (30/40 minutes)
- Fast results (4/12 sessions).

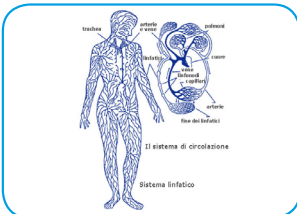

Elimination of adiposity through the circulatory and lymphatic system.



Developed by the father of ultrasounds.



# OSMOLIPO

## INVERSE ULTRASONIC OSMOSIS

**ACTIONS:** • FAT REDUCTION



3 Transducers to treat different depth of fat :


- Flat (2/3cm)
- Concave (3/4cm)
- Convex (1/2 cm)
- Mini



Mini handprobe to treat small areas. Ideal for double chin.



Effective therapeutic results also for the reduction of lipomas (medical procedures)

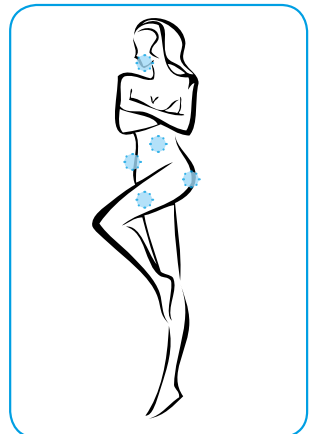

Treatable areas:

- Buttocks • Thighs • Culottes de chevalle
- Knees • Arms • Abdomen • Double chin

MADE IN ITALY


# OSMOLIPO

## INVERSE ULTRASONIC OSMOSIS

### THE TECHNOLOGY:

ITS Group research and development department created OSMOLIPO, a device to eliminate localized fat. OSMOLIPO is the result of extensive research between private industry and university laboratories in the field of non-invasive technologies for the treatment of imperfections of the body: the most efficient and technologically advanced device available today in medicine and aesthetic. The evolution in the use of ultrasound on the human body has allowed the development of 'inverse ultrasonic osmosis', an innovative system of ultrasonic waves that allows the progressive increase of the permeability of the cytoplasmic membrane causing leakage of the contents (triglycerides and fatty acids) with consequent elimination of adipocytes through the lymphatic and circulatory system.

The results are visible in few sessions (4-12) and it is possible to reshape and remodel any area of the body.

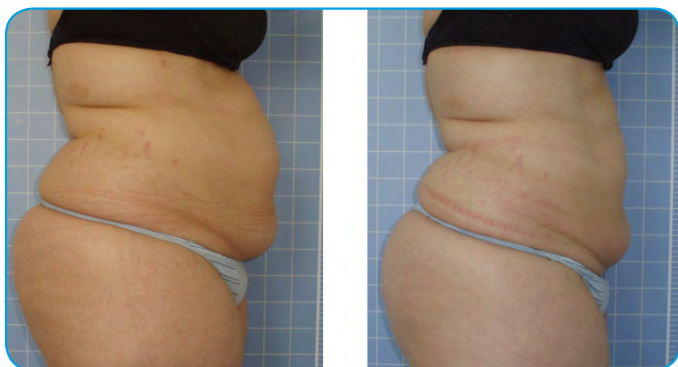

### CONTRAINDICATIONS:



- Diabetes
- Renal failure / dialysis
- Prescribed drugs
- Obesity
- Anti-inflammatory (chemical drugs)
- Cholesterol problems
- Metal implants over the treatment area
- Menstrual period (only for abdominal treatment)
- Pregnancy
- Pacemaker wearers

### INFORMATION:



- Number of sessions: 4/12
- Session interval : weekly
- Session duration : 30/40 minutes

### TECHNICAL SPECIFICATIONS:

Model:	OSMOLIPO
Supply :	110-240V 50-60Hz
Power consumption:	300W
Frequency:	33 KHz
Power consumption:	3W/cm2
Classification 93/42/CEE:	II A
Classification IEC 60601-1:	Classe I tipo B
Display:	Touch Screen
Size:	L 440 mm, P 360 mm, H 250 mm
Weight:	12 Kg.
Compliance and Standards:	Temperature environment: 15°C - 35°C, moisture 30% - 60%, pressure 86kPa - 106kPa
Storage and transport:	Temperature environment: 5°C - 45°C, moisture 30% - 70%, pressure 86kPa - 106kPa



Via Del Cerchio 1 - 40012 Calderara di Reno (BO) - Italy  
Tel. (039) 051 6467153 - Fax: (039) 051 0822612  
info@itsgroupsrl.com - www.itsgroupsrl.com



ISO 9001  
ISO 13485



Rev. 16.02